




Pressure Mapping Analysis

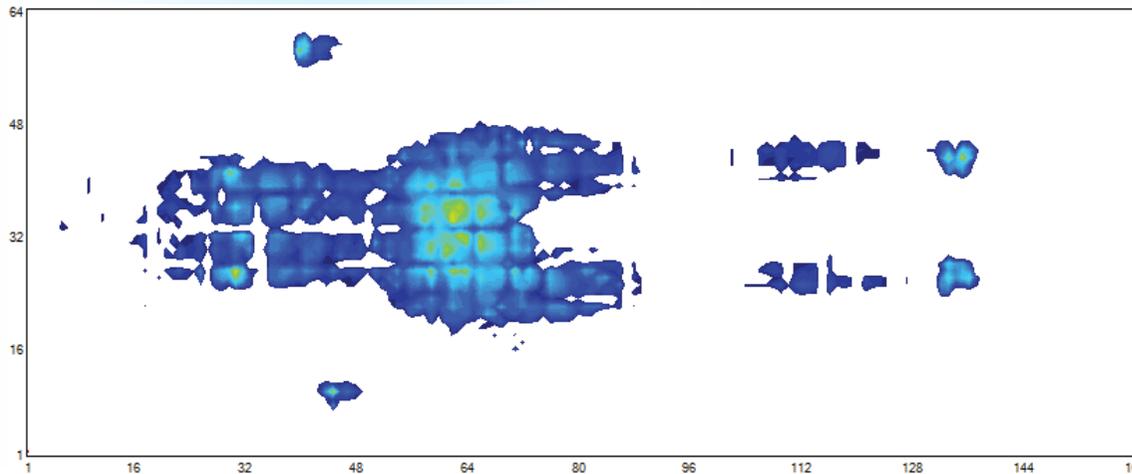

Female 65 kg Hospital bed laid flat
Avg. Press. 18.49
Peak Press. 41.52
Min. Press. 10.34
Area (cm²) 2412.9

mmHg	
56.12	0
54.54	0
52.97	0
51.39	0
49.81	0
48.23	0
46.65	0
45.07	0
43.49	0
41.91	0
40.34	3
38.76	5
37.18	12
35.6	20
34.02	14
32.44	16
30.86	22
29.29	36
27.71	34
26.13	43
24.55	57
22.97	80
21.39	83
19.81	89
18.24	121
16.66	117
15.08	149
13.5	194
11.92	193
10.34	210
0.0	8742
10.34 - 56.12	
Sensor Group	Value
*X100:64.160:02 S0123	
Avg Pres.	18.49
Peak Pres.	41.52
Min Pres.	10.34
Area (cm ²)	2412.9

RMS

Replacement Foam Mattress

Description

Foam mattress, castellated on both sides to aid in air circulation when in use, and to help the mattress contour to the bed frame when profiled. This mattress has a 4 way turning cycle. Two way stretch wipe clean vapour permeable P.U. cover and zip on three sides.

Dimensions

(L)1980 x (W)910 x (H)150 mm

Warranty

3 years

Max User Weight

190.5Kg/30 stones

Fire Retardancy

BS 7177 : 2008
 1996 Medium Hazard Ignition sources 0.1.5. (Crib 5), BS 6807 : 1990